

From discovery to application: Enabling scalable, reproducible targeted proteomics through standardized Evosep workflows

Djordje Vasiljević, Iulia Macavei, Magnus Huusfeldt, Stine R. Rasmussen, Maurine Fucito, Stoyan Stoychev, Nicolai Bache, Dorte Bekker-Jensen

Evosep Biosystems, Denmark

Highlights

- Evosep Proteomics enables a high-performance Sample to Separation solution for standardized and scalable mass spectrometry.
- Robust ecosystem for scalable proteomics with RT reproducibility below 1 s RT SD and peptide CVs below 5% using the Evosep Eno 200 SPD method.
- Highly robust assay sensitivity from Evokits with reproducible quantification from sub-nanogram to microgram protein lysate input levels.


From discovery to targeted application

From discovery to application through a standardized ecosystem

Despite the large number of identified biomarker candidates, only a limited subset has been translated into robust assays for clinical research. Targeted LC–MS enables this transition by converting discovery-derived targets into quantitative verification and validation workflows. Achieving this at scale requires standardized and controlled processes that ensure consistent execution across samples and studies.


Evosep Proteomics addresses this by integrating sample preparation and LC–MS into a unified ecosystem. Evokits, Evosep Lupo, and Evosep Eno form a seamless workflow, with Evtip Pure linking preparation to analysis.

Together, these components operate as a single controlled system, enabling consistent transfer from discovery datasets to scalable targeted workflows.

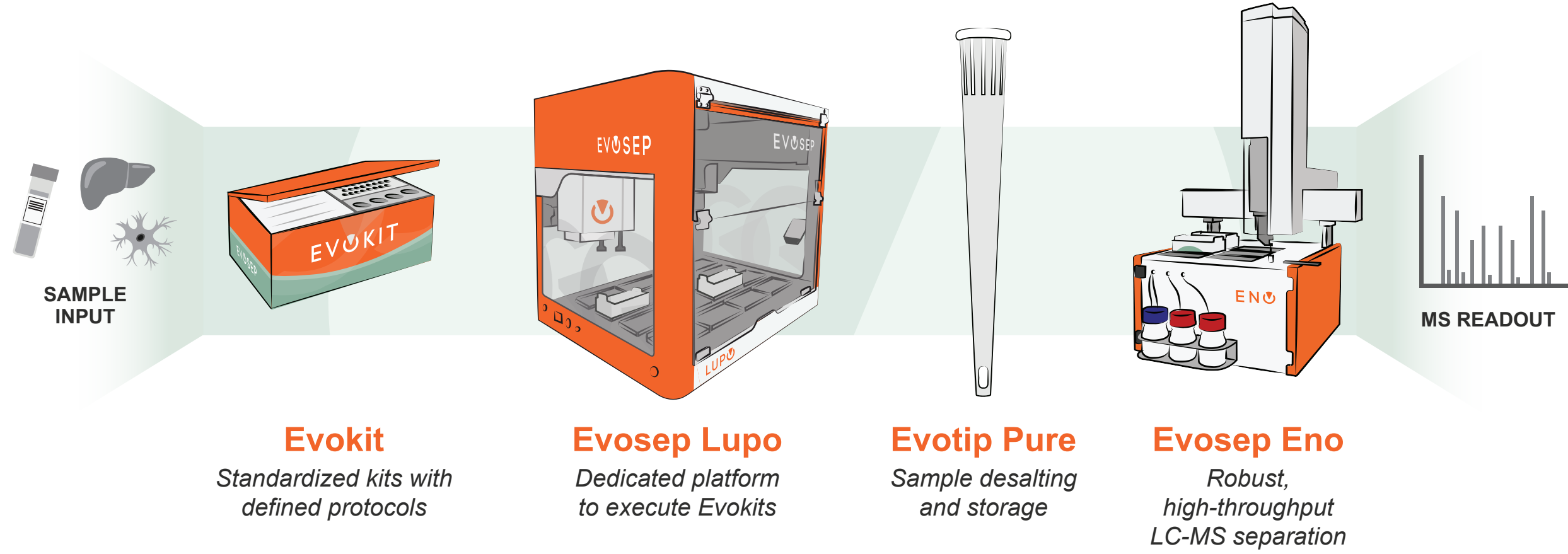

Built-in controls and walk-away automation for robust results


The workflow enables automated processing from protein extract to purified digest loaded on Evtip using Evokits, which include all required reagents, consumables and Evtips. Execution on Evosep Lupo incorporates built-in traceability and quality controls, while predefined workflows and software interface minimize user intervention. Integrated QC samples monitor both sample preparation and LC–MS performance throughout the study. Defined peptide panels track digestion efficiency, workflow recovery, and instrument stability, ensuring controlled and reproducible targeted analysis.

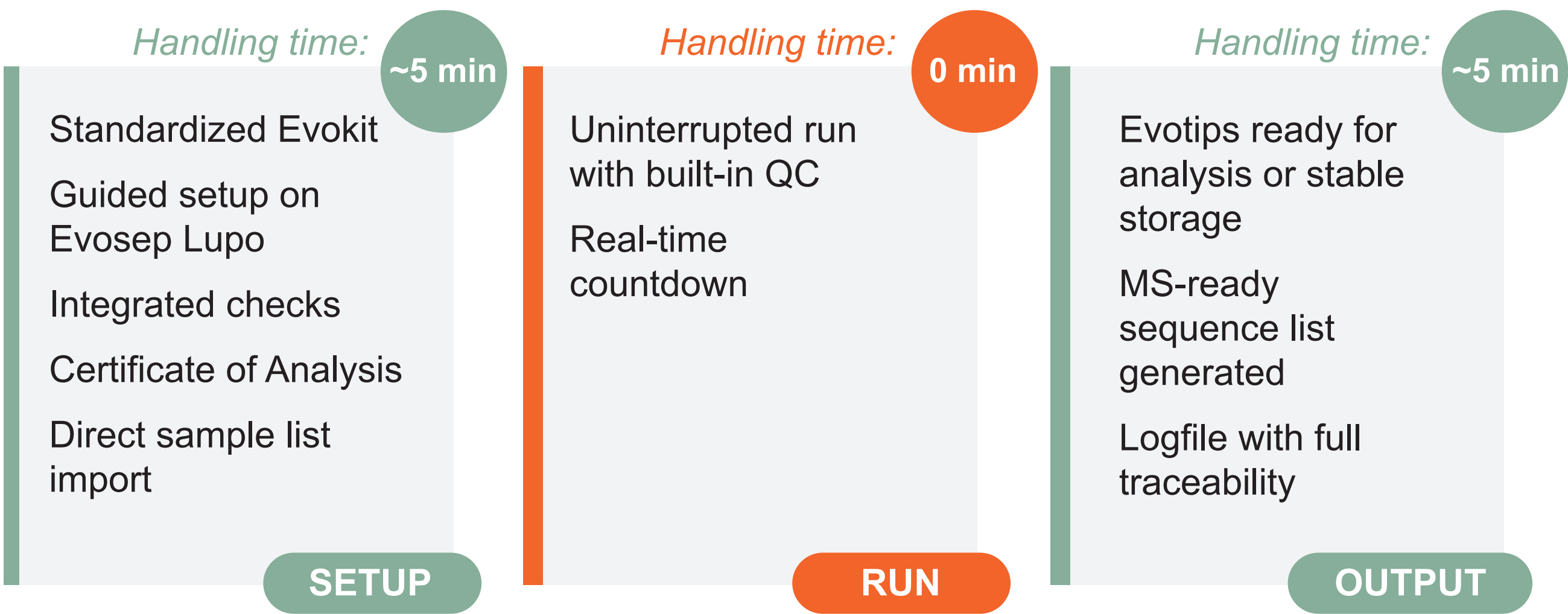

Controlled, consistent results across every workflow step

High sensitivity for targeted workflows


The sensitivity of the Evokit Digest was measured in MRM mode on an Evosep Eno 200 SPD coupled to a Waters Xevo TQ-Absolute XR. Extracted ion chromatograms and calibration curves show high linear signal response across ~3 orders of magnitude ($R^2 \approx 0.996\text{--}0.999$) using 1 to 1024 ng protein lysate with the Evokit, with 30% of the digest loaded on the Evtip for analysis. Notably, great quantification is obtained with very little starting material, demonstrating robust workflow sensitivity. The recovery remains within ~80–120% across the range, confirming accurate and reproducible quantification from sub-nanogram to microgram levels.

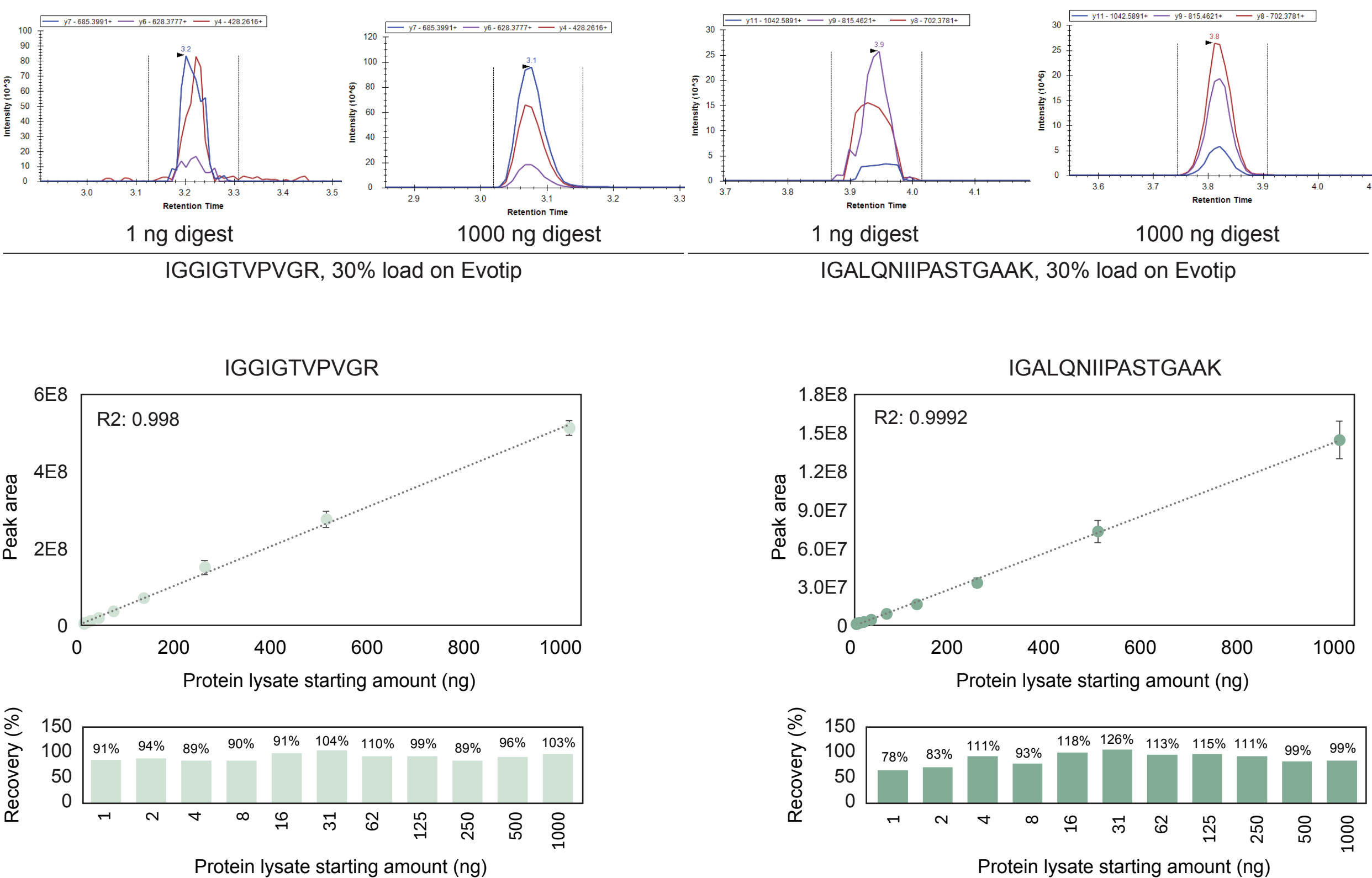

Sensitivity of Evokit Digest is shown with two different peptides, measured on a Evosep Eno 200 SPD, coupled to a Waters Xevo TQ-Absolute XR (n=8). Extracted ion chromatograms (top), calibration curves (middle) and recovery (bottom).

Seamless transition from discovery to targeted analysis

Evosep Proteomics enables a seamless transition from discovery to targeted analysis. This is facilitated by standardized sample preparation with Evokits using the Evosep Lupo, reproducibly generating consistent peptide features across studies as well as LC-MS analyses. Three independent digests prepared on different days were analyzed using two discovery MS systems and one QQQ MS, each coupled to a different Evosep Eno instrument. Across all conditions, the same peptide features were consistently detected with high precision (CV <10%) and retention time stability (RSD), demonstrating robustness independent of sample preparation day, LC setup, and MS platform.

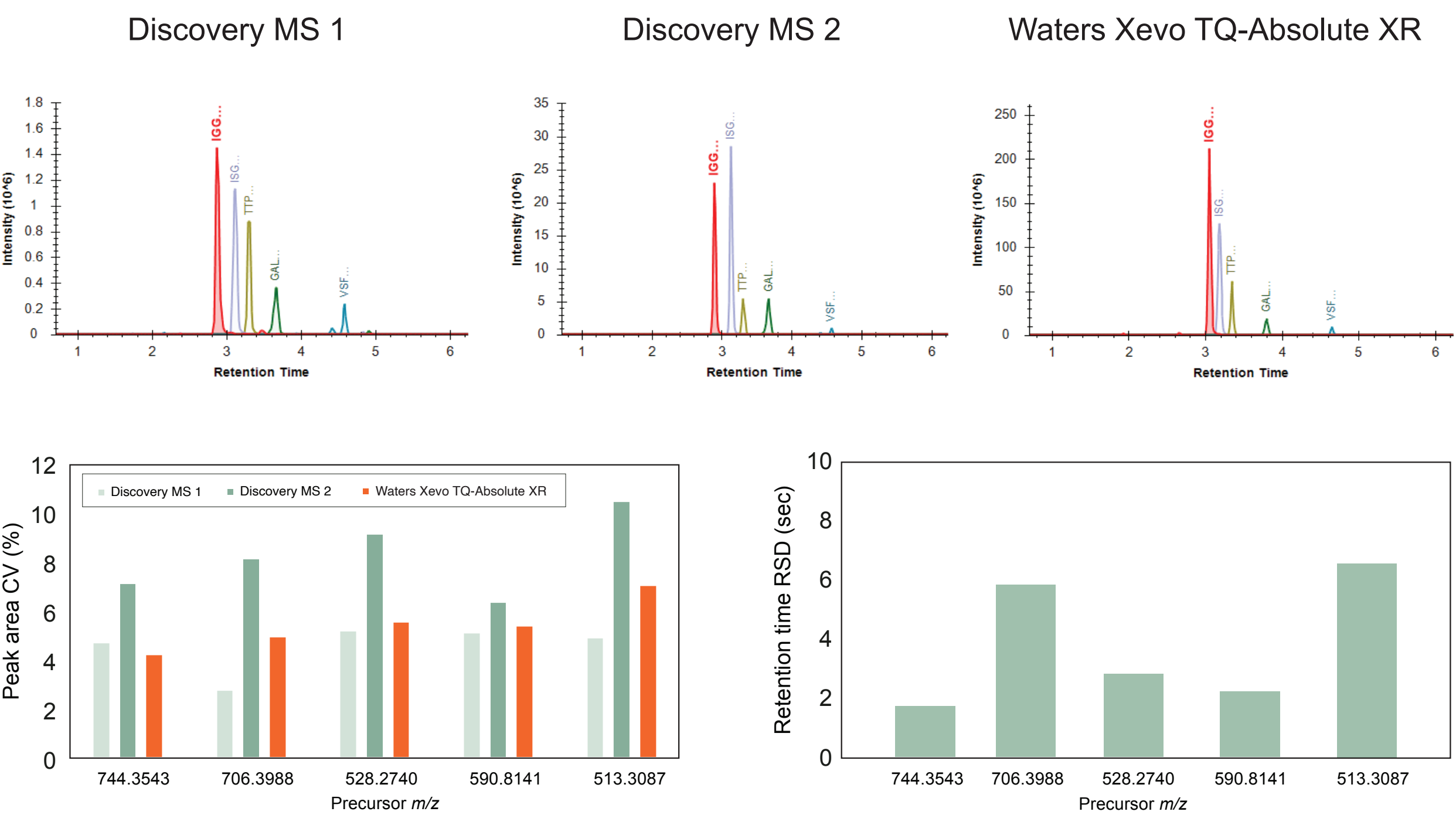

Extracted ion chromatograms, workflow precision and retention time reproducibility across the Discovery and Targeted MS platforms on different days, Evosep Eno 200 SPD. 1 µg K562 lysate digested with Evokit. 30% peptide loaded on Evtip (n=6).

Sample preparation at scale with Evokits and Evosep Lupo

Built-in quality controls incorporated with each Evokit

Each Evokit Digest includes both process control samples and system suitability controls, which are used to monitor the sample preparation and to verify LC–MS performance over time. This setup enables continuous quality assessment of pre-analytical handling and instrument stability throughout the analytical process. While, the QC samples are processed in the last column of each plate, they are analyzed across all samples during LC-MS. Predefined run order setups are easily integrated through an easy-to-use software interface.

